

daisygrip the reusable tourniquet


intuitive, self-finding
magnetic buckle



buckle and band
able to be completely
disinfected

innovative, wipeable
silicone band for patient
comfort and safety

award winning

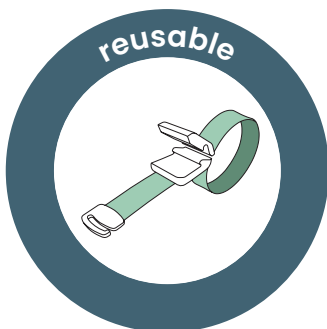

why daisygrip?



daisygrip is a reusable tourniquet, intended to slow the flow of blood for the purpose of intravenous access. As a **Red Dot Design Award Winner***, the daisygrip has taken an essential medical device, and redeveloped it in line with the real-world need for infection control, sustainability and economy.

daisygrip unique features:

- Construction materials mean the entire device can be wiped for effective disinfection immediately after use.
- Band made of skin-friendly materials, which also prevent pinching.
- In vitro studies and publications demonstrate that the device both picks up less contamination during use and is able to be more effectively disinfected than conventional tourniquets.⁴
- Intuitive, self-finding closure can be operated with a single hand.

daisygrip – a better tourniquet

An Australian study⁵ found that – compared to a fabric tourniquet – the daisygrip was considered more environmentally sustainable, able to be cleaned more effectively and easier to clean between patients.

It was the preferred choice by users. The facility also estimated significant annual cost savings, as well as a reduction in landfill and carbon output compared to disposables used in some areas.

Contact us for a copy of the poster!⁶



daisygrip provides:

- A better experience for patient and user
- Higher standard of infection control
- Costs savings compared to conventional methods
- Increased sustainability practices

current evidence has highlighted that:

- More than **70%** of tourniquets exhibit contamination.¹
- Construction materials for both conventional reusable and single-use tourniquets have been shown to **pick up and transfer microorganisms**.²
- A hospital in New Zealand found various levels of **contamination** on conventional tourniquets, with the highest levels on those located on phlebotomy trolleys after the ward rounds. Crucially, the facility disinfects all tourniquets overnight.³

references:

*Red Dot Design Award: daisygrip (red-dot.org)

¹ "Health professionals' practices related with tourniquet use during peripheral venipuncture: a scoping review" (de Sousa Salgueiro-Oliveira et al. Rev. Latino-Am. Enfermagem) <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC6528630/pdf/0104-1169-rlae-27-e3125.pdf>

² "Methicillin resistant Staphylococcus aureus contamination of phlebotomy tourniquets and faucets" (Abeywickrama et al. Ceylon Medical Journal) <https://cmj.sjoiinfo/articles/10.4038/cmj.v63i1.8627>

³ "Quantifying patient bacterial exposure risk from reusable phlebotomy tourniquets in a NZ secondary level hospital" (Schauer and Hammer. Journal of Inf. Prevention) https://www.ncbi.nlm.nih.gov/pmc/articles/PMC5074164/pdf/10.1177_1757177415600242.pdf

⁴ "Reduced bacterial contamination rates detected on silicone tourniquets compared to conventional tourniquets" (Graham et al. BMC Infectious Diseases) Reduced bacterial contamination rates detected on silicone tourniquets compared to conventional tourniquets in clinical routine | BMC Infectious Diseases | Full Text (biomedcentral.com)

⁵ <https://acipconference.com.au/comparison-of-single-patient-use-and-multi-patient-use-fabric-tourniquets-with-a-silicone-reusable-tourniquet-daisygrip-in-an-australian-tertiary-hospital/>

⁶ <https://acipc2023-c75494.epresentor.com.au/posters?tid=9700&preview=aHR0cHM6Ly9hY2lwZWwMjMTYzc-INDk0LmVwcmVzW50ZiXiuY29tLmFjL2NsaWVudHMvNzlvODMwL3NIYmJpc3Npb25zLzE3NjA4My9wb3N0ZXlvcG-9zdGVyX3dlYjE2OTg4OTIzNDUvaW5kZXguaHRtbA==&theme=Going%20Green%20%26%20Sustainability&posterId=176083&author=Kate%20Ryan&indiDownload=0&tv=148>

Distributed by:

Tristel™

Tristel Solutions Limited
Unit 1b, Lynx Business Park, Fordham Road,
Snailwell, Cambridgeshire, CB8 7NY, UK
T +44 (0) 1638 721500 – E mail@tristel.com
W www.tristel.com

UK Person Responsible
RIVINGTON SOLICITORS LIMITED
Equitable House, Suite 207, 7 General
Gordon Sq., London, SE18 6FH, UK
Sushil Gaikwad: rivington@daisygrip.de



Refer to the
User Guide
for full product
instructions.

**Scan to
find out
more**

